



# Arena 2P

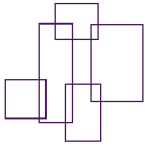

## Suggested fitting range:

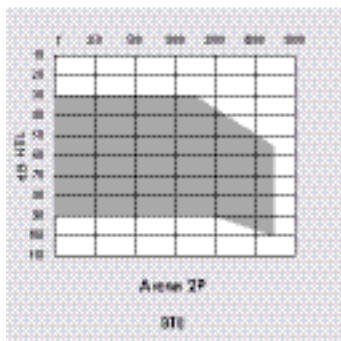

## Arena 2P (Power) BTE

### Technical Description:

The Arena 2P is a power, 100% digital signal processing BTE instrument. This two channel instrument offers features such as a low battery warning and Microphone Noise Reduction for improved listening and wearability. The Arena 2P BTE includes three standard trimmers: Low Cut (NH), AGC-I (TK), and AGC-O.

### Fitting Recommendation:

The Arena 2P BTE is an excellent solution for moderate to severe hearing losses.

### Standard Features:

- Fully digital WDRC amplifier
- Three adjustable trimpot controls (NH, TK, and AGC-O)
- Two channels
- Adjustable broadband AGC-O
- Microphone noise reduction
- Volume control
- Telecoil
- M-T-O switch
- Low battery alert tone
- Direct audio input
- Compatible with FM, MLX, and other wireless technology
- Beige, grey, or brown housing

### Available Options:

- Transparent eyeglass adapter
- Audio shoe
- Small earhook with damper
- Volume control cover

ANSI S3.22-1996

MAXIMUM POWER OUTPUT	Peak	dB OSPL90	130
	HF-average	dB OSPL90	125
GAIN (Input at 50 dB)	Peak	dB	64
	HF-average	dB	57
REFERENCE TEST GAIN		dB	48
FREQUENCY RESPONSE		Hz	100-5600
COMPRESSION	Attack time	ms	VARIABLE
	Release time	sec	VARIABLE
TOTAL HARMONIC DISTORTION	Maximum	%	<10
EQUIVALENT NOISE LEVEL	Maximum	dB SPL	<33
INDUCTION COIL SENSITIVITY	HFA-SPLITS@31.6 mA/m		
	Left/Right	dB SPL	107/109
BATTERY DRAIN (maximum)	at RTG., Input 65 dB SPL	mA	0.9
TYPICAL BATTERY LIFE	Type 13 zinc air	hrs	331

5010 Cheshire Ln. N.  
Suite 2  
Plymouth, MN 55446  
800-876-1141  
Fax: 763-553-9129  
www.rexton-online.com

**REXTON**

# Arena 2P (Power) BTE

