


## Targa PRO FREE



### Suggested fitting range:

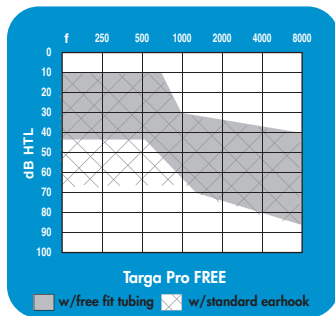

## Targa PRO FREE BTE Spec Sheet

### Technical Description:

The Targa PRO FREE BTE is an intermediate level, 100% DSP circuit with two channels, two adjustable AGC-I compression channels, and three independently programmable memories. Targa's unique technology allows for extremely quick and accurate fittings. This BTE features a compact case and the thin FREE fit tubing caters to cosmetically-conscious consumers. The Targa PRO FREE Open Fitting System offers the hearing health professional and the consumer a choice between open tubing and a traditional fit with an attached earmold. The Targa PRO FREE BTE also offers directional microphones, an enhanced feedback cancellation system, adjustable signal tone loudness, power-on-delay and low battery signal. Targa products are programmed using the Connex and RexFit fitting software.

### Fitting Recommendation:

The Targa PRO FREE BTE is an excellent choice for those with mild to moderate sloping or high frequency hearing losses as well as for those patients who prefer a non-occluding ear mold.

### Standard Features:

- 2 channels
- 3 Memories
- Directional Microphone
- Feedback Cancellation
- Telecoil
- 2 Adjustable AGC-I
- 2 Slope Controls
- Crossover Frequency Adjustment
- Power-On-Delay
- Low Battery Warning
- Push Button On/Off
- 13A Battery
- Beige, grey and brown housing

### Available Options:

- Eyeglass adapter
- Small Earhook

ANSI S3.22-2003			STANDARD EARHOOK	Targa FREE FIT TUBING
MAXIMUM POWER OUTPUT	Peak	dB OSPL90	124	123
	HF-average	dB OSPL90	122	115
GAIN (Input at 50 dB)	Peak	dB	55	47
	HF-average	dB	51	41
REFERENCE TEST GAIN		dB	45	38
FREQUENCY RESPONSE		Hz	100-6300	100-5800
COMPRESSION (AGC-O @ -21 DB)	Attack time	ms	VARIABLE	VARIABLE
	Release time	ms	VARIABLE	VARIABLE
TOTAL HARMONIC DISTORTION	Maximum	%	7	4
EQUIVALENT NOISE LEVEL		dB SPL	19	19
INDUCTION COIL SENSITIVITY	HFA-SPLITS@31.6 mA/m			
	Left/Right	dB SPL	104/102	97/94
BATTERY DRAIN	at RTG., Input 65 dB SPL	mA	0.8	0.8
BATTERY LIFE	Type 13 zinc air	hrs	~280	~280

5010 Cheshire Ln. N.  
Suite 2  
Plymouth, MN 55446  
800-876-1141  
Fax: 763-553-9129  
[www.rexton-online.com](http://www.rexton-online.com)

**REXTON**

# Targa PRO FREE BTE Spec Sheet

