


Calibra 1 PRO P



Suggested fitting range:

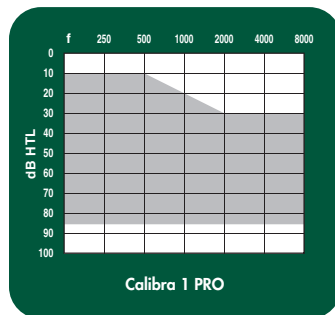

Calibra 1 PRO P BTE

Technical Description:

The Calibra 1 PRO P (power) directional BTE is a premier-level 100% DSP circuit with four channels, four independently programmable memories with one memory dedicated solely for audio input. This BTE features an endless volume control wheel and Automatic Microphone Management. Other outstanding features of the Calibra 1 PRO P include DataLogging, Feedback Cancellation, Equalized Directional Microphone, and four personalized choices for Environment and Noise Management; speech, speech in noise, noise and music. Rexton's easy-to-use software helps you program this superior technology by incorporating new and improved fitting protocols in a familiar programming platform.

Fitting Recommendation:

The Calibra 1 PRO Power BTE is an excellent choice for mild to severe hearing loss of many configurations including flat, sloping, high frequency, and reverse curve.

Standard Features:

- 4 Channels
- 4 Memories, 1 dedicated for audio input
- Volume Control
- Automatic Directional Microphone
- Feedback Cancellation
- Personalized Noise Reduction Management
- DagaLogging
- 4 AGC Adjustments
- Power On Delay
- Memory Selection Tones
- Push Button On/Off
- Telecoil
- Compatible with wireless products, cell-phones, FM systems, and MLX
- Beige, grey, brown and metallic champagne housing color
- 13A battery type

Available Options:

- Eyeglass adapter
- Small Earhook
- Audio Shoe/Integrated Battery Door

ANSI S3.22-2003

MAXIMUM POWER OUTPUT	Peak	dB OSPL90	130
	HF-average	dB OSPL90	126
GAIN (Input at 50 dB)	Peak	dB	65
	HF-average	dB	58
REFERENCE TEST GAIN		dB	48
FREQUENCY RESPONSE		Hz	100-6600
COMPRESSION	Attack time	ms	VARIABLE
	Release time	ms	VARIABLE
TOTAL HARMONIC DISTORTION	Maximum	%	7
EQUIVALENT NOISE LEVEL		dB SPL	16
INDUCTION COIL SENSITIVITY	HFA-SPLITS@31.6 mA/m		
	Left/Right	dB SPL	105/109
BATTERY DRAIN	at RTG., Input 65 dB SPL	mA	0.9
BATTERY LIFE	Type 13 zinc air	hrs	~250

Rexton, Inc.

5010 Cheshire Ln. N., Suite 2
Plymouth, MN 55446

800-876-1141

Fax: 763-553-9129

www.rexton-online.com

REXTON

Calibra 1 PRO P BTE

