


## Targa PLUS P



### Suggested fitting range:

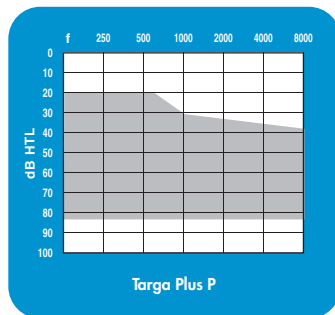

## Targa PLUS P (power) BTE Spec Sheet

### Technical Description:

The Targa PLUS P BTE instrument is an intermediate level, 100% DSP circuit with four channels, four adjustable compression channels and four independently programmable memories with one memory dedicated solely for audio input. The Targa PLUS P BTE features a directional microphone, adaptive noise reduction, output limiting and user friendly controls such as power-on-delay, low battery warning, programmable on/off switch, and adjustable signal tones.

### Fitting Recommendation:

The Targa PLUS P BTE is an excellent choice for a variety of hearing losses, including moderate to severe losses of many configurations including flat, sloping, high frequency and reverse curve.

### Standard Features:

- 4 Channels
- 4 AGC-I Adjustments
- Directional Microphones
- Adaptive Noise Reduction
- Volume Control
- Telecoil
- Adjustable Slope Controls
- Adjustable Crossover Frequency
- Feedback Cancellation
- Microphone Noise Reduction
- 4 Programmable Memories with 1 memory dedicated to Audio Input
- Power-On-Delay
- Low Battery Warning
- Adjustable Signal Tone Loudness
- Beige, grey, brown and metallic champagne housing
- 13A Battery Size

### Optional Features:

- Eyeglass Adapter
- Small Earhook
- Audio Shoe/Integrated Battery Door

### ANSI S3.22-2003

MAXIMUM POWER OUTPUT	Peak	dB OSPL90	132
	HF-average	dB OSPL90	126
GAIN (Input at 50 dB)	Peak	dB	65
	HF-average	dB	57
REFERENCE TEST GAIN		dB	49
FREQUENCY RESPONSE		Hz	100-6400
COMPRESSION	Attack time	ms	VARIABLE
	Release time	ms	VARIABLE
TOTAL HARMONIC DISTORTION	Maximum	%	7
EQUIVALENT NOISE LEVEL		dB SPL	16
INDUCTION COIL SENSITIVITY	HFA-SPLITS@31.6 mA/m		106/110
	Left/Right	dB SPL	
BATTERY DRAIN	at RTG., Input 65 dB SPL	mA	0.8
BATTERY LIFE	Type 13 zinc air	hrs	~280

### Rexton, Inc.

5010 Cheshire Ln. N., Suite 2  
 Plymouth, MN 55446  
 800-876-1141  
 Fax: 763-553-9129  
[www.rexton-online.com](http://www.rexton-online.com)

# REXTON

# Targa PLUS P (power) BTE Spec Sheet

